






JDENTALCARE

just smile



MADE IN ITALY

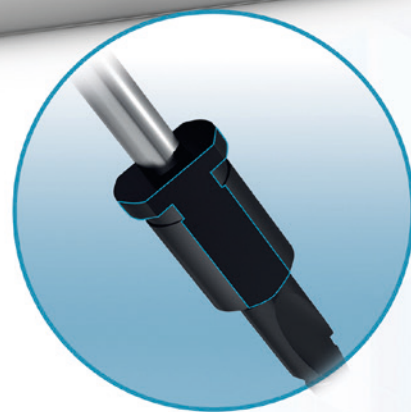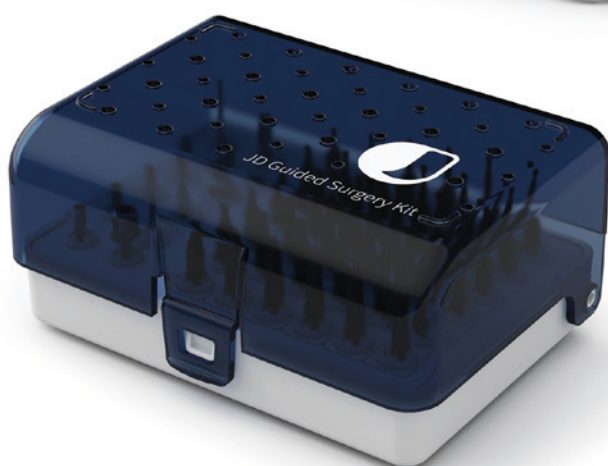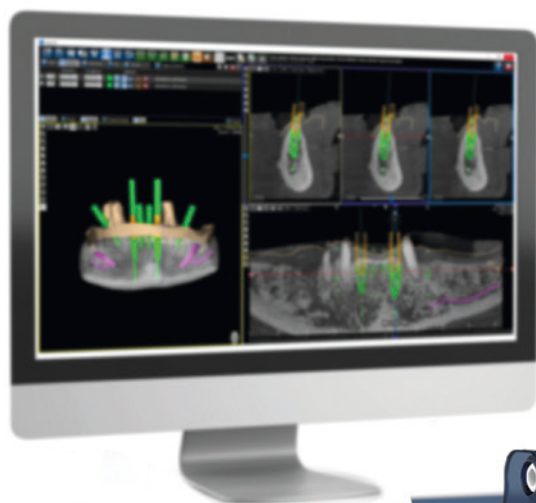

JD Guided Surgery Kit

Friend of your hand

LET YOUR HAND BE GUIDED BY THE JD GUIDED SURGERY KIT

The ideal solution to perform a minimally invasive surgery with more precision, greater predictability of outcomes and with few simple guided surgery steps.

NEW DESIGN FOR BETTER SURGERY

All JD Guided Drills have a non cylindrical design with two lateral cuts which allows:



- Lower friction in the sleeves so as to avoid overheating
- Better external irrigation during the implant site preparation



PRECISION

GUIDED IMPLANT DEPTH CONTROL

High precision in surgical procedures: implants are placed to the planned depth using a surgical guide and the guided surgery instruments are used in combination with a single sleeve of Ø 5.01mm for a secure implant depth control.

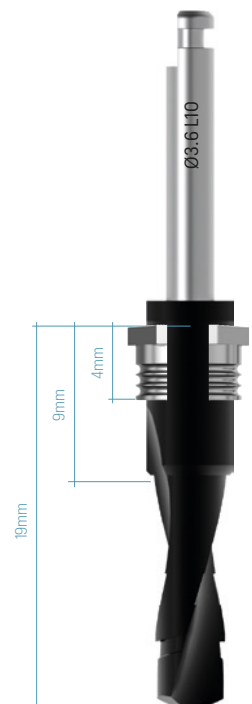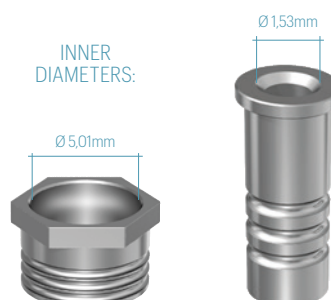


Note: all measurements in mm

SIMPLICITY

ONE KIT FOR ALL JD IMPLANT LINES

Simplified insertion with just one kit for all JD implant lines from 2.75 to Ø 5.0mm and from lengths 6 to 15mm.

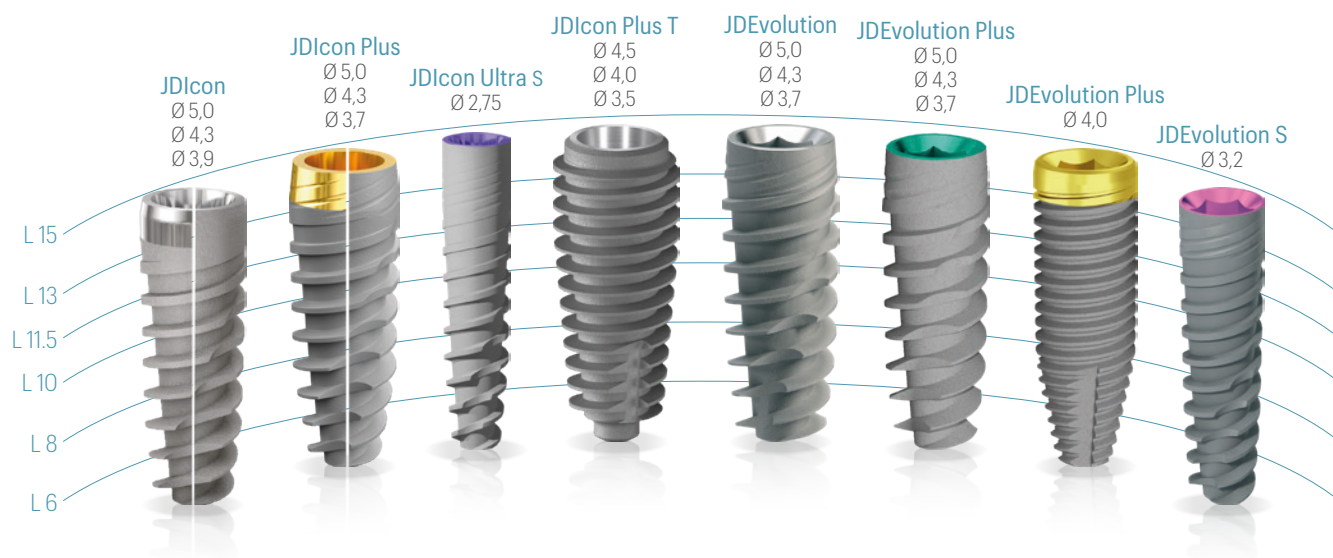


COMPATIBILITY

USER-FRIENDLY SOFTWARE DESIGN

A single kit compatible with the main Software programs for guided surgery present on the market.

PRODUCTS

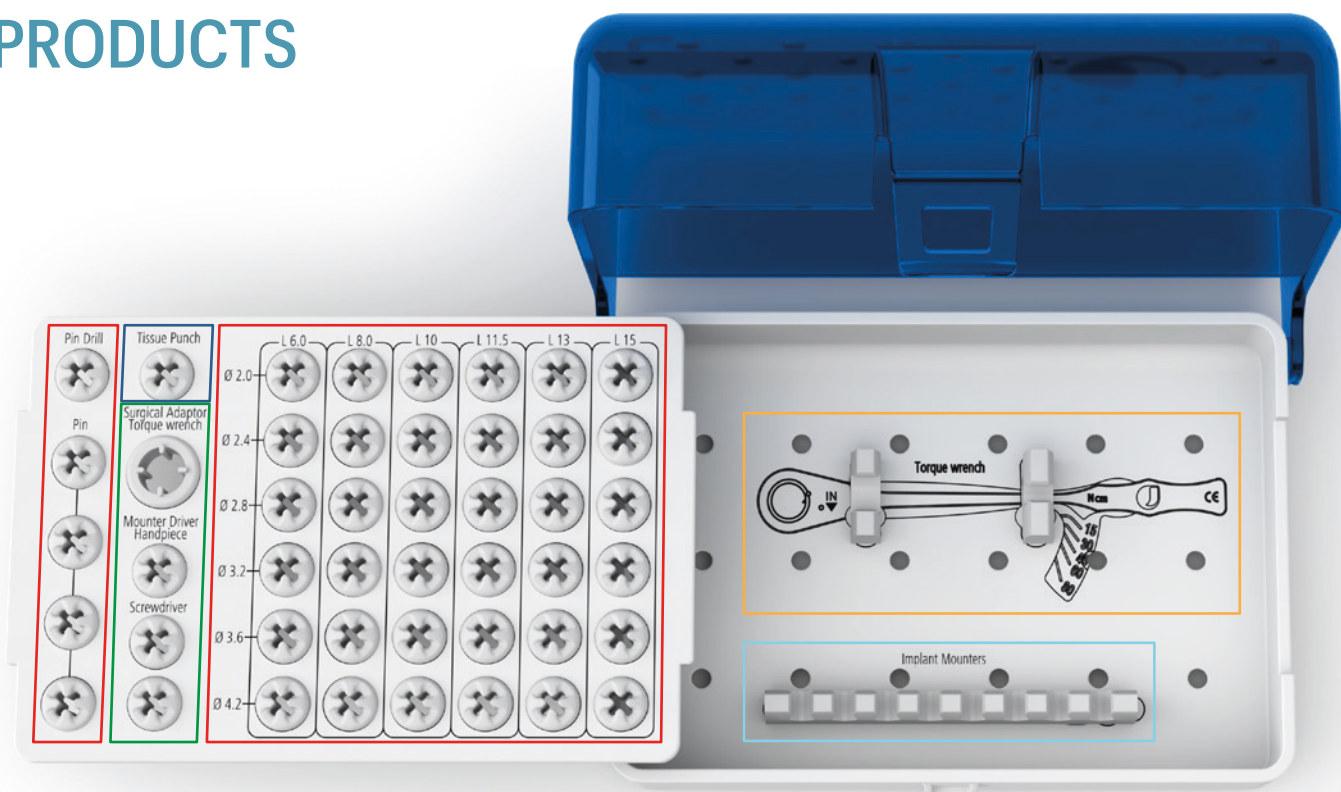

JD Torque



JDTW
Torque Wrench JDTorque

Implant and Prosthetic Drivers



EVIDG
Driver Body for Guided Surgery



JDID116*
Implant Driver for Guided Surgery JDEvolution Plus



JDID111*
Implant Driver for Guided Surgery JDEvolution S



EVSDP25
Prosthetic Driver L 25



EVSDPF25
Prosthetic Driver for JDTorque L 25



JDTWA
Surgical Adaptor for JDTorque

Tissue Punch



JDTP35
Tissue Punch Ø 3.5

Guided implant Mounters



EVGMC
Implant Mounter for JDEvolution



ESGMC
Implant Mounter for JDEvolution S



ICGMC.
Implant Mounter for JDIcon Plus



ISGMC
Implant Mounter for JDIcon Ultra S



EVGMLC
Implant Mounter Long for JDEvolution



ESGMLC
Implant Mounter Long for JDEvolution S



ICGMLC.
Implant Mounter Long for JDIcon Plus



ICGMC
Implant Mounter for JDIcon

Guided Drills



JDGD20-060
Guided Drill
Ø 2.0 L 6.0



JDGD20-080
Guided Drill
Ø 2.0 L 8.0



JDGD20-100
Guided Drill
Ø 2.0 L 10.0



JDGD20-115
Guided Drill
Ø 2.0 L 11.5



JDGD20-130
Guided Drill
Ø 2.0 L 13.0



JDGD20-150
Guided Drill
Ø 2.0 L 15.0

Note: all measurements in mm

*To be ordered separately

Guided Drills



JDGD24-060
Guided Drill
Ø 2.4 L 6.0



JDGD24-080
Guided Drill
Ø 2.4 L 8.0



JDGD24-100
Guided Drill
Ø 2.4 L 10.0



JDGD24-115
Guided Drill
Ø 2.4 L 11.5



JDGD24-130
Guided Drill
Ø 2.4 L 13.0



JDGD24-150
Guided Drill
Ø 2.4 L 15.0



JDGD28-060
Guided Drill
Ø 2.8 L 6.0



JDGD28-080
Guided Drill
Ø 2.8 L 8.0



JDGD28-100
Guided Drill
Ø 2.8 L 10.0



JDGD28-115
Guided Drill
Ø 2.8 L 11.5



JDGD28-130
Guided Drill
Ø 2.8 L 13.0



JDGD28-150
Guided Drill
Ø 2.8 L 15.0



JDGD32-060
Guided Drill
Ø 3.2 L 6.0



JDGD32-080
Guided Drill
Ø 3.2 L 8.0



JDGD32-100
Guided Drill
Ø 3.2 L 10.0



JDGD32-115
Guided Drill
Ø 3.2 L 11.5



JDGD32-130
Guided Drill
Ø 3.2 L 13.0



JDGD32-150
Guided Drill
Ø 3.2 L 15.0



JDGD36-060
Guided Drill
Ø 3.6 L 6.0



JDGD36-080
Guided Drill
Ø 3.6 L 8.0



JDGD36-100
Guided Drill
Ø 3.6 L 10.0



JDGD36-115
Guided Drill
Ø 3.6 L 11.5



JDGD36-130
Guided Drill
Ø 3.6 L 13.0



JDGD36-150
Guided Drill
Ø 3.6 L 15.0



JDGD42-060
Guided Drill
Ø 4.2 L 6.0



JDGD42-080
Guided Drill
Ø 4.2 L 8.0



JDGD42-100
Guided Drill
Ø 4.2 L 10.0



JDGD42-115
Guided Drill
Ø 4.2 L 11.5



JDGD42-130
Guided Drill
Ø 4.2 L 13.0



JDGD42-150
Guided Drill
Ø 4.2 L 15.0



JDGD100*
Guided Drill
Countersink



JDGDP
Guided Surgery
Precision Drill



JDPIN
Fixation Pin
(x4)

Sleeves



JDBGB
Template Sleeve for
Guided Surgery Peek



JDBG
Template Sleeve
for Guided Surgery



JDBGP
Pin Sleeve

Note: all measurements in mm

*To be ordered separately

DRILLING PROTOCOLS

The following drilling protocols are suggested to be used for the insertion of JDentalCare dental implants in case of Guided Surgery technique, using the original instruments and surgical guides manufactured by JDentalCare. Deviations in respect to the following procedures are possible depending on specific clinical situations.

JDEvolution, JDEvolution Plus, JDIcon, JDIcon Plus, JDIcon Plus T, JDOcta					
		Site preparation in maxilla		Site preparation in mandible	
IMPLANT DIAMETER	IMPLANT LENGTH	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
Ø 3,7 / 3,9	L6	/	/	/	/
	L8	2,4 L6 - L8 2,8 L6 3,2 L6	2,4 L6 - L8 2,8 L6 3,2 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 3,6 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 - L8 3,6 L6
	L10	2,4 L6 - L8 - L10 2,8 L6 3,2 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 - L10 3,6 L6
	L11,5	2,4 L6 - L8 - L10 - L11,5 2,8 L6 3,2 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 - L11,5 3,6 L6
	L13	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 - L13 3,6 L6
	L15	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6
	L18	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6
Ø 4,0 / 4,3	L6	2,4 L6 2,8 L6 3,2 L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6 4,2 L6
	L8	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 3,6 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 3,6 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 3,6 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 - L8 3,6 L6 - L8 4,2 L6
	L10	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 - L10 3,6 L6 - L8 - L10 4,2 L6
	L11,5	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 3,6 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 - L11,5 3,6 L6 - L8 - L10 - L11,5 4,2 L6
	L13	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 3,6 L6 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 - L13 3,6 L6 - L8 - L10 - L11,5 - L13 4,2 L6
	L15	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 3,6 L6 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 - L11,5 - L13 - L15 4,2 L6
	L18	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 3,6 L6 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 3,6 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 - L11,5 - L13 - L15 4,2 L6

Note: All measurements in mm

<i>JDEvolution, JDEvolution Plus, JDIcon, JDIcon Plus, JDIcon Plus T, JDOcta</i>					
		Site preparation in maxilla		Site preparation in mandible	
IMPLANT DIAMETER	IMPLANT LENGTH	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
Ø 5,0	L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6 4,2 L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6 4,2 L6	2,4 L6 2,8 L6 3,2 L6 3,6 L6 4,2 L6	/
	L8	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 - L8 3,6 L6 4,2 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 - L8 3,6 L6 4,2 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6 - L8 3,6 L6 - L8 4,2 L6	/
	L10	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 - L10 3,6 L6 4,2 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 - L10 3,6 L6 - L8 4,2 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6 - L8 - L10 3,6 L6 4,2 L6	/
	L11,5	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 - L11,5 3,6 L6 4,2 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 - L11,5 3,6 L6 - L8 - L10 4,2 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6 - L8 - L10 - L11,5 3,6 L6 - L8 4,2 L6	/
	L13	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 - L13 3,6 L6 - L8 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 - L13 3,6 L6 - L8 - L10 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 - L11,5 - L13 3,6 L6 - L8 - L10 4,2 L6	/
	L15	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 4,2 L6	/
	L18	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 - L11,5 4,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6 - L8 - L10 - L11,5 - L13 - L15 3,6 L6 - L8 - L10 - L11,5 4,2 L6	/

<i>JDIcon Ultra S</i>				
IMPLANT DIAMETER	IMPLANT LENGTH	SOFT BONE TYPE IV	MEDIUM BONE TYPE II-III	DENSE BONE TYPE I
Ø 2,75	L6	/	/	/
	L8	2,0 L6 - L8	2,0 L6 - L8 2,4 L6 - L8	2,0 L6 - L8 2,4 L6 - L8 2,8 L6
	L10	2,0 L6 - L8 - L10	2,0 L6 - L8 - L10 2,4 L6 - L8 - L10	2,0 L6 - L8 - L10 2,4 L6 - L8 - L10 2,8 L6
	L11,5	2,0 L6 - L8 - L10 - L11,5	2,0 L6 - L8 - L10 - L11,5 2,4 L6 - L8 - L10 - L11,5	2,0 L6 - L8 - L10 - L11,5 2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8
	L13	2,0 L6 - L8 - L10 - L11,5 - L13	2,0 L6 - L8 - L10 - L11,5 - L13 2,4 L6 - L8 - L10 - L11,5 - L13	2,0 L6 - L8 - L10 - L11,5 - L13 2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8

<i>JDEvolution S</i>					
		Site preparation in maxilla		Site preparation in mandible	
IMPLANT DIAMETER	IMPLANT LENGTH	SOFT BONE	MEDIUM-DENSE BONE	SOFT BONE	MEDIUM-DENSE BONE
Ø 3,2	L6	/	/	/	/
	L8	2,4 L6 - L8 2,8 L6	2,4 L6 - L8 2,8 L6 3,2 L6	2,4 L6 - L8 2,8 L6 3,2 L6	2,4 L6 - L8 2,8 L6 - L8 3,2 L6
	L10	2,4 L6 - L8 - L10 2,8 L6	2,4 L6 - L8 - L10 2,8 L6 3,2 L6	2,4 L6 - L8 - L10 2,8 L6 3,2 L6	2,4 L6 - L8 - L10 2,8 L6 - L8 - L10 3,2 L6
	L11,5	2,4 L6 - L8 - L10 - L11,5 2,8 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 3,2 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 2,8 L6 - L8 - L10 - L11,5 3,2 L6
	L13	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6
	L15	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 3,2 L6	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 2,8 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6

JD EXTRA DRILLS KIT

Note: For implants lengths from 18mm to 26mm use the JD Extra Drills Kit.
The kit contains:

JDGD20-180	Guided Drill Ø 2.0 L 18.0
JDGD24-180	Guided Drill Ø 2.4 L 18.0
JDGD20-200	Guided Drill Ø 2.0 L 20.0
JDGD24-200	Guided Drill Ø 2.4 L 20.0
JDGD20-220	Guided Drill Ø 2.0 L 22.0
JDGD24-220	Guided Drill Ø 2.4 L 22.0
JDGD20-240	Guided Drill Ø 2.0 L 24.0
JDGD24-240	Guided Drill Ø 2.4 L 24.0
JDGD20-260	Guided Drill Ø 2.0 L 26.0
JDGD24-260	Guided Drill Ø 2.4 L 26.0

Note: all measurements in mm



DRILLING PROTOCOLS

JD <i>P</i> terygo		
IMPLANT DIAMETER	IMPLANT LENGTH	
Ø 3,3	L13	2,0 L6 - L8 - L10 - L11,5 - L13 2,4 L6 2,8 L6 (only for dense bone)
	L15	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 2,4 L6 2,8 L6 (only for dense bone)
	L18	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 2,4 L6 2,8 L6 (only for dense bone)
	L20	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 2,4 L6 2,8 L6 (only for dense bone)
Ø 4,0	L13	2,4 L6 - L8 - L10 - L11,5 - L13 3,2 L6
	L15	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 3,2 L6
	L18	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 3,2 L6
	L20	2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 3,2 L6

Note: All measurements in mm

JD <i>N</i> asal		
IMPLANT DIAMETER	IMPLANT LENGTH	
Ø 4,0	L20	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6
	L22	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6
	L24	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 - L24 2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 - L24 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6
	L26	2,0 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 - L24 - L26 2,4 L6 - L8 - L10 - L11,5 - L13 - L15 - L18 - L20 - L22 - L24 - L26 2,8 L6 - L8 - L10 - L11,5 - L13 3,2 L6 - L8 - L10 3,6 L6